


# Circumferential RF Solution

Only 360° heating coil RF catheter for  
incompetent perforator and tributary  
veins\* that avoids the need  
for repeated quadrant ablations\*\*

**Venclose Maven™**

Perforator Catheter

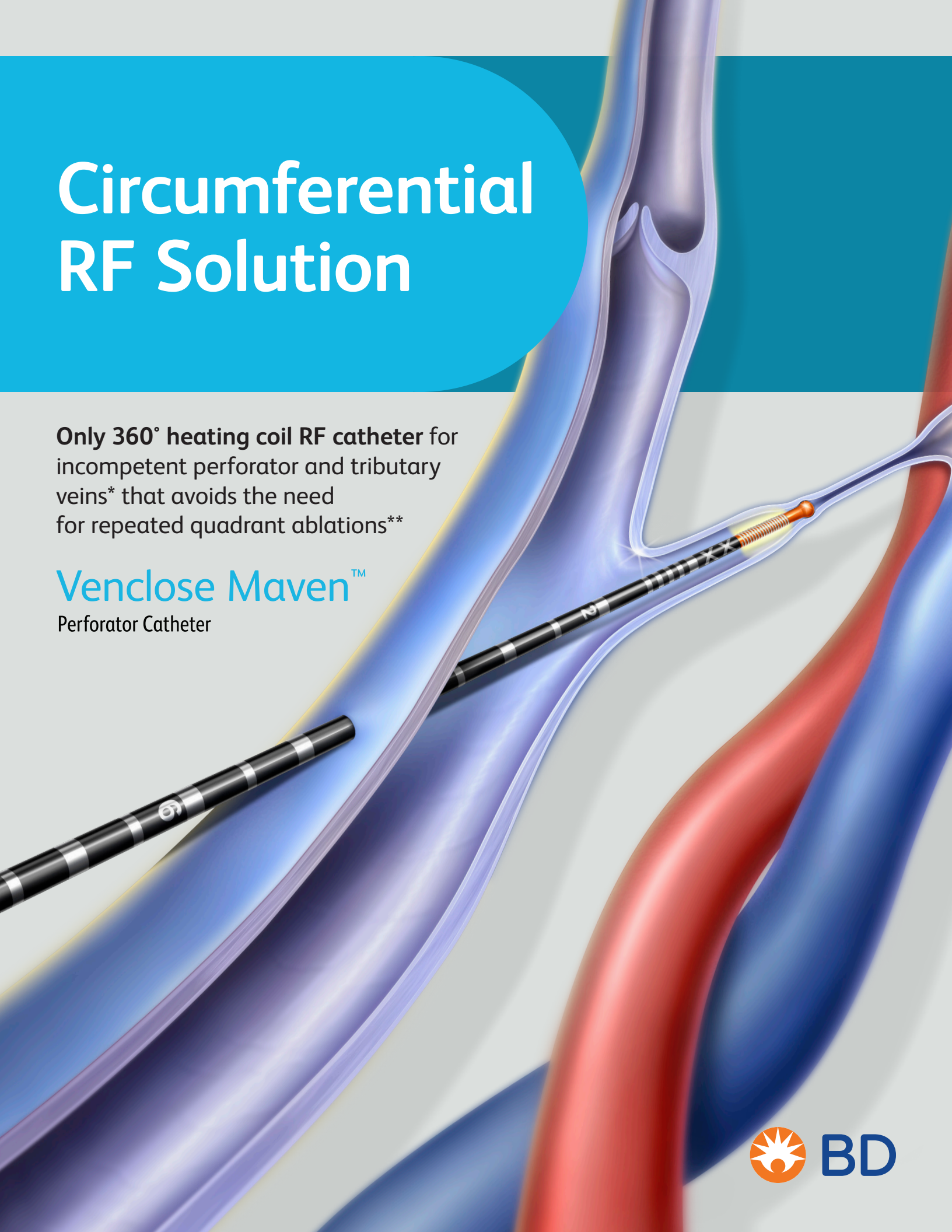

# A Modernized, 360° Approach



**LOWERS TOTAL  
ABLATION TIME**  
COMPARED TO CLOSURERFS™ STYLET\*\*

**HELPS TO MINIMIZE  
INVASIVENESS**  
WITH LOW-PROFILE 6F CATHETER

**AVOIDS THE NEED  
FOR REPEATED  
QUADRANT ABLATIONS\*\***



**360°**

CIRCUMFERENTIAL RESISTIVE  
HEATING FOR 20-SECONDS PER  
TREATMENT CYCLE




**FLEXIBLE**

40 CM FLEXIBLE CATHETER DESIGN  
SUPPORTS NAVIGATION FOR VARYING  
VEIN LENGTHS AND ANATOMIES



**CURVED**

CATHETER TO HELP  
FACILITATE POSITIONING

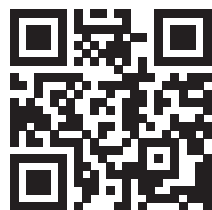

**Learn More  
By Visiting  
Venclose.com**

**Venclose Maven™**  
Perforator Catheter

\*As of September 2024

\*\*The Venclose Maven™ Perforator Catheter uses 0.5 cm segmental ablation with a 20-second ablation cycle. The ClosureRFS™ Stylet uses pull back or focal ablation. Both ClosureRFS™ Stylet ablation methods require an initial quadrant ablation treatment time of 4 minutes. The pull back method requires additional quadrant ablation treatments for 1 minute per 1 cm, and the focal method requires additional quadrant ablation of subsequently treated adjacent vein segments for 4 minutes each. The ClosureRFS™ Stylet information was obtained from the manufacturer's procedure guidelines.

The Venclose Maven™ Catheter is intended for endovascular coagulation of blood vessels in patients with perforator and tributary vein reflux. The Venclose Maven™ Catheter is contraindicated in patients with thrombus in the vein segment to be treated. Potential adverse events include, but are not limited to the following: vessel perforation, skin discoloration, nerve injury, temporary paresthesia, thrombosis, deep vein thrombosis, phlebitis, hematoma, infection, skin burn, pulmonary embolism, and pain. **Please consult product labels and instructions for use for indications, contraindications, hazards, warnings, and precautions.**

## Ordering Information

### Perforator Catheter

Description	Product Codes
Venclose Maven™ Perforator Catheter (40 cm)	<input type="checkbox"/> VC056F

### Generator & Accessories

Description	Product Codes
Venclose™ RF Generator	<input type="checkbox"/> VCRFG1
Venclose™ Procedure Pack (No Access)	<input type="checkbox"/> VCPK
12G Angiocath™ IV Catheter	<input type="checkbox"/> 382277
7 cm Micro Introducer Sheath	<input type="checkbox"/> NIS02A
Venclose™ System Foot Pedal	<input type="checkbox"/> VCFP1
Venclose™ System US Power Cord	<input type="checkbox"/> VCPCB

**bd.com**

BD, Tempe, AZ, USA, 1 800 321 4254

BD, the BD Logo, Angiocath, Venclose, and Venclose Maven are trademarks of Becton, Dickinson and Company or its affiliates. © 2024 BD. All Rights Reserved. © 2024 Illustration by Mike Austin. BD-54911v2

